

## Product Change Notification

Subject: Add degree markings to the polarizer dial on the BL-ZW1 Backlight Pad

Products Affected:

BL-ZW1 Dino-Lite Backlight Pad

Description of Change:

The polarizer dial of the BL-ZW1 now have laser inscribed markings (fig.1) at 5 degree increments for referencing the degree of polarization. When the degree indicator on the BL-ZW1 is aligned (fig.2) with the magnification indicator on the Dino-Lite, 0 and 180 degree positions on the BL-ZW1 indicates there is no polarization, of which the polarizer of the BL-ZW1 is in parallel to the polarizer at the sensor side of the Dino-Lite.

Note: The polarizer of BL-ZW1 will take effect only when it is paired with a Dino-Lite model with polarizer.

Effective Date

Immediately

ANMO

For any questions regarding this notification, please contact your AnMo contact person.

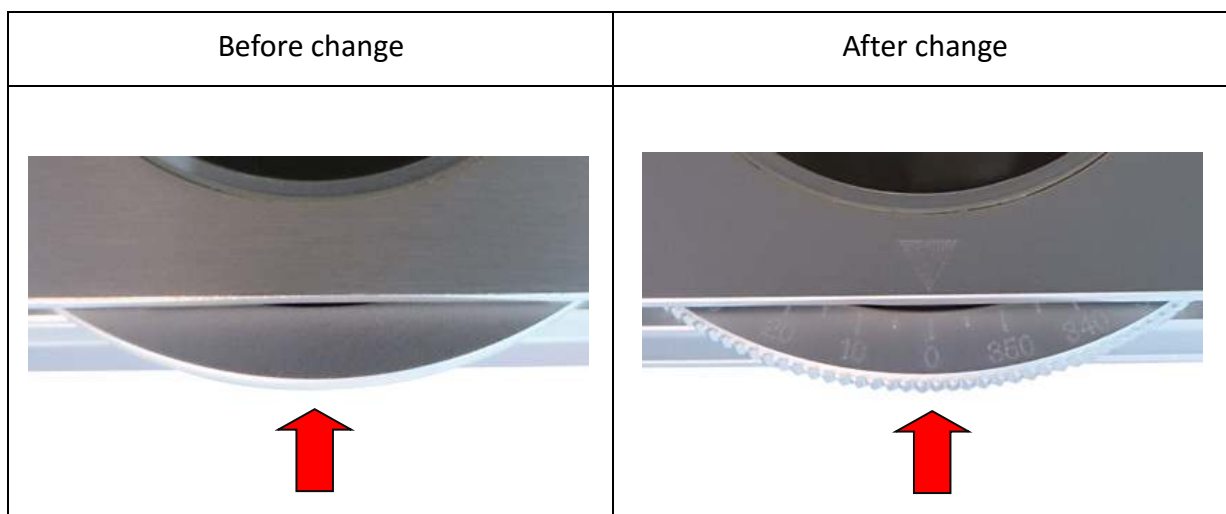


Fig.1



**Fig.2.** Align the degree indicator of the BL-ZW1 to the magnification indicator of a Dino-Lite with polarizer.